

Technical Data

Transmission ratio	2.8 : 1
Maximum speed	2.700 – 4.200 rpm
Cannula stroke	2.0 mm
Weight, complete	490 g

Ref

Conform Cannula Handpiece, complete	5077nou
Handle complete for Liposuction	5107
Silicone sealing CHP (5 pcs)	40378

+ NOUVAG AG

St. Gallerstrasse 25
9403 Goldach
Switzerland
T+41 (0) 71 846 66 00
info@nouvag.com

+ NOUVAG GmbH

Schulthaißstrasse 15
78462 Konstanz
Germany
T+49 (0) 7531 1290-0
F+49 (0) 7531 1290-12
info-de@nouvag.com

nouvag.com

NOUVAG+

A Masterpiece
even better


THE NEW LIPOSUCTION HANDPIECE **MARK II**

©NOUVAG AG • 67173 • V20230818 • All rights reserved

THE NEW LIPOSUCTION HANDPIECE **MARK II**

Manufactured in Switzerland

Optimizations for **MARK II**

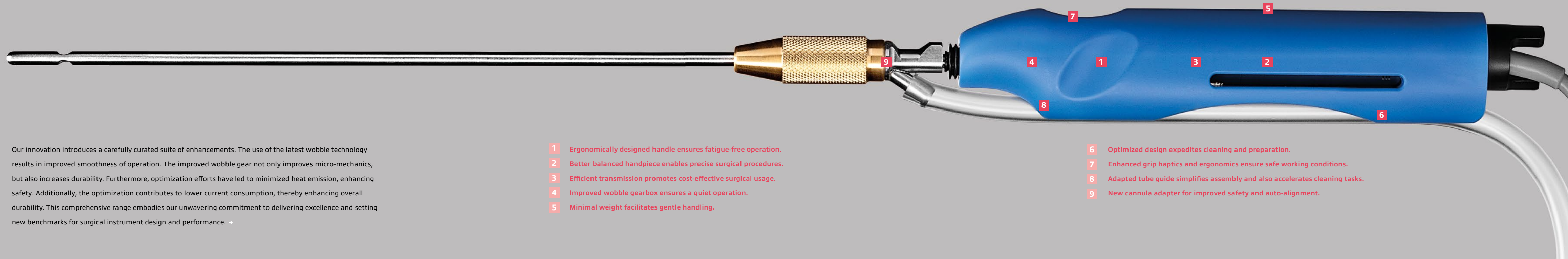

Our innovation introduces a carefully curated suite of enhancements. The use of the latest wobble technology results in improved smoothness of operation. The improved wobble gear not only improves micro-mechanics, but also increases durability. Furthermore, optimization efforts have led to minimized heat emission, enhancing safety. Additionally, the optimization contributes to lower current consumption, thereby enhancing overall durability. This comprehensive range embodies our unwavering commitment to delivering excellence and setting new benchmarks for surgical instrument design and performance. →

- 1** Ergonomically designed handle ensures fatigue-free operation.
- 2** Better balanced handpiece enables precise surgical procedures.
- 3** Efficient transmission promotes cost-effective surgical usage.
- 4** Improved wobble gearbox ensures a quiet operation.
- 5** Minimal weight facilitates gentle handling.

- 6** Optimized design expedites cleaning and preparation.
- 7** Enhanced grip haptics and ergonomics ensure safe working conditions.
- 8** Adapted tube guide simplifies assembly and also accelerates cleaning tasks.
- 9** New cannula adapter for improved safety and auto-alignment.